

**This content is not available in your preferred language.**

The content is shown in another available language. Your browser may include features that can help translate the text.

# Import CSV Spreadsheet Files into LabVIEW

Updated Mar 26, 2019

Reported In 

## Software

- LabVIEW Base
- LabVIEW Full
- LabVIEW Professional

## Other

- Microsoft Office
- Microsoft Excel


## Issue Details

I have a Microsoft Excel spreadsheet that I would like to import into LabVIEW, but I do not want to use ActiveX calls to read the data. Is there another way to read in this type of file?


## Solution

One way to read an Excel spreadsheet into LabVIEW is by saving it as a `.csv` file, which is a feature built into Excel. To save a spreadsheet in Excel as a `.csv` and import it into LabVIEW, follow the below steps:

1. In Excel, select **Save As** and in
  - Microsoft Office 2010 and older: Select **Other Formats** (shown below)
  - Microsoft Office 2013: select a location to save.

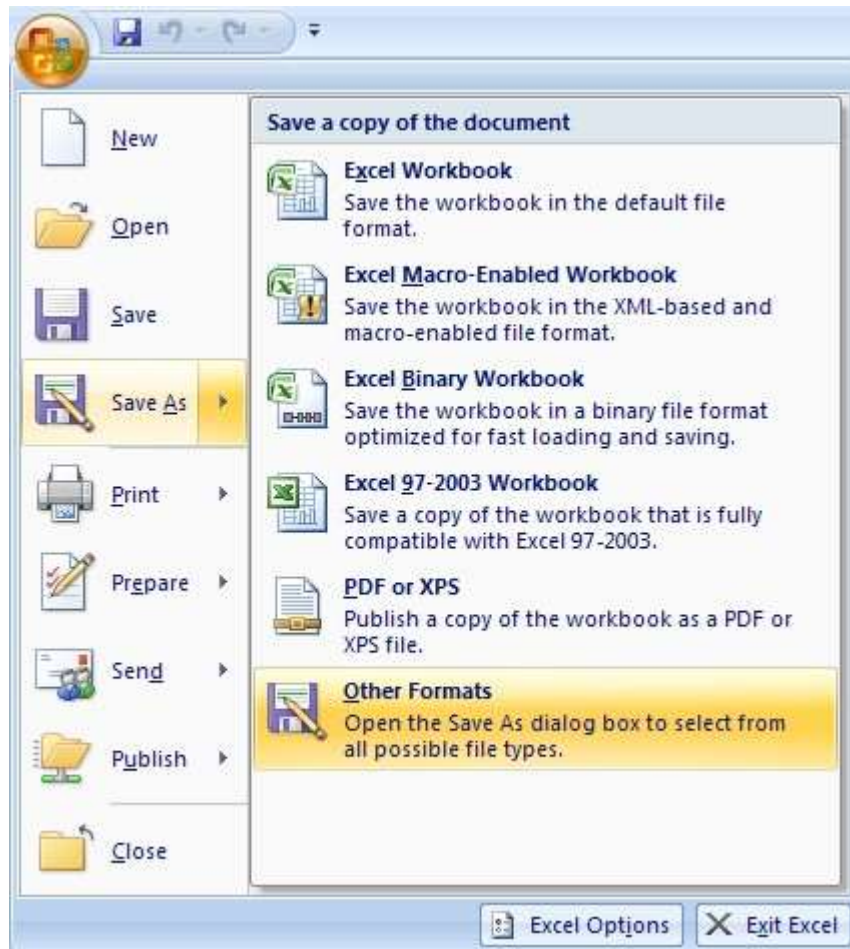


2. In the **Save As** type drop down menu, select **CSV (Comma delimited)(\*.csv)**.

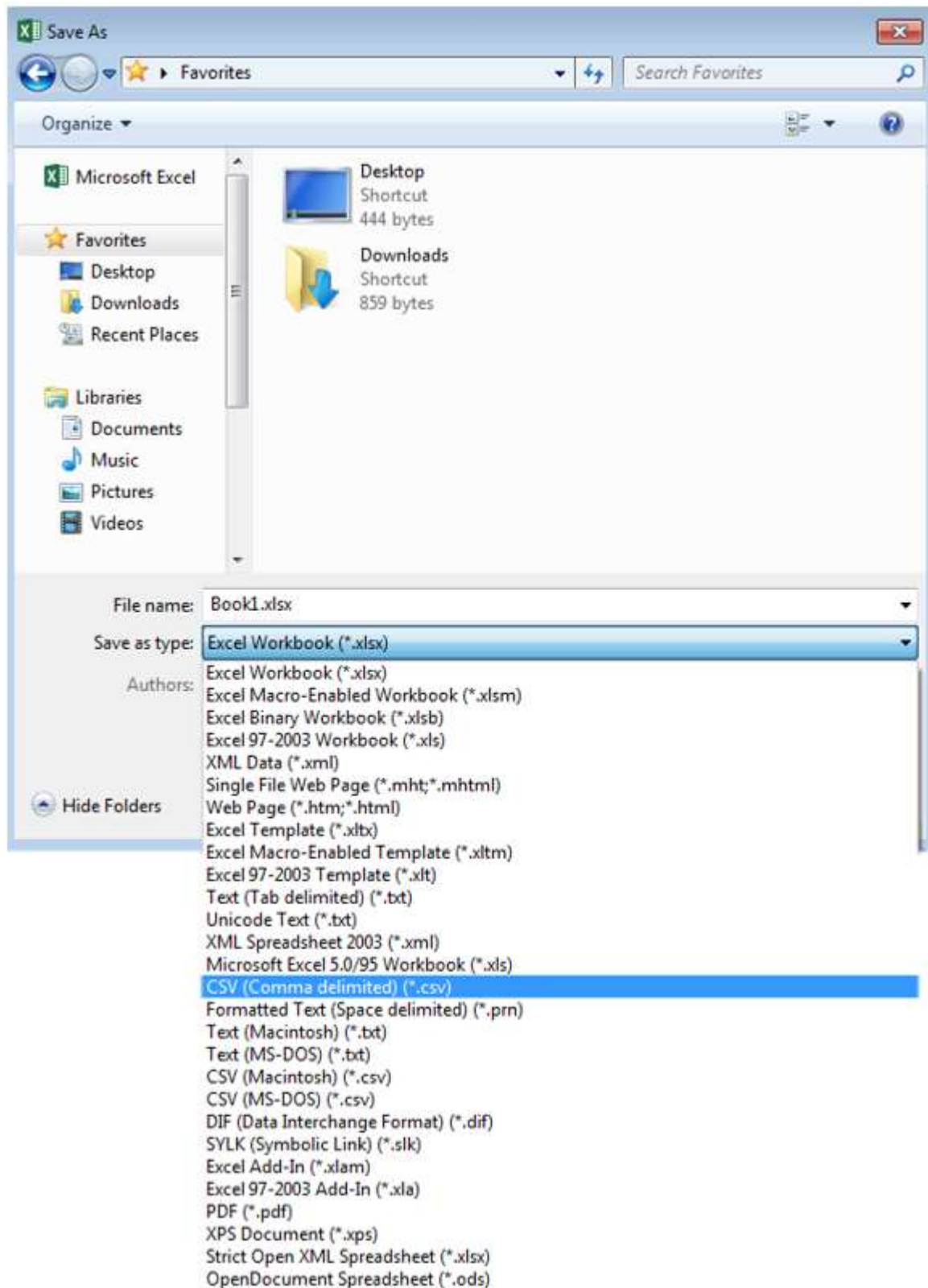


**Note:** Two caveats--using the .csv file format does not support multiple sheets, and some other Excel features may not be compatible. Excel will display two warnings about these limitations when you save the file.

3. In LabVIEW, bring up the Functions Palette and select **File I/O » Read Delimited Spreadsheet.vi**

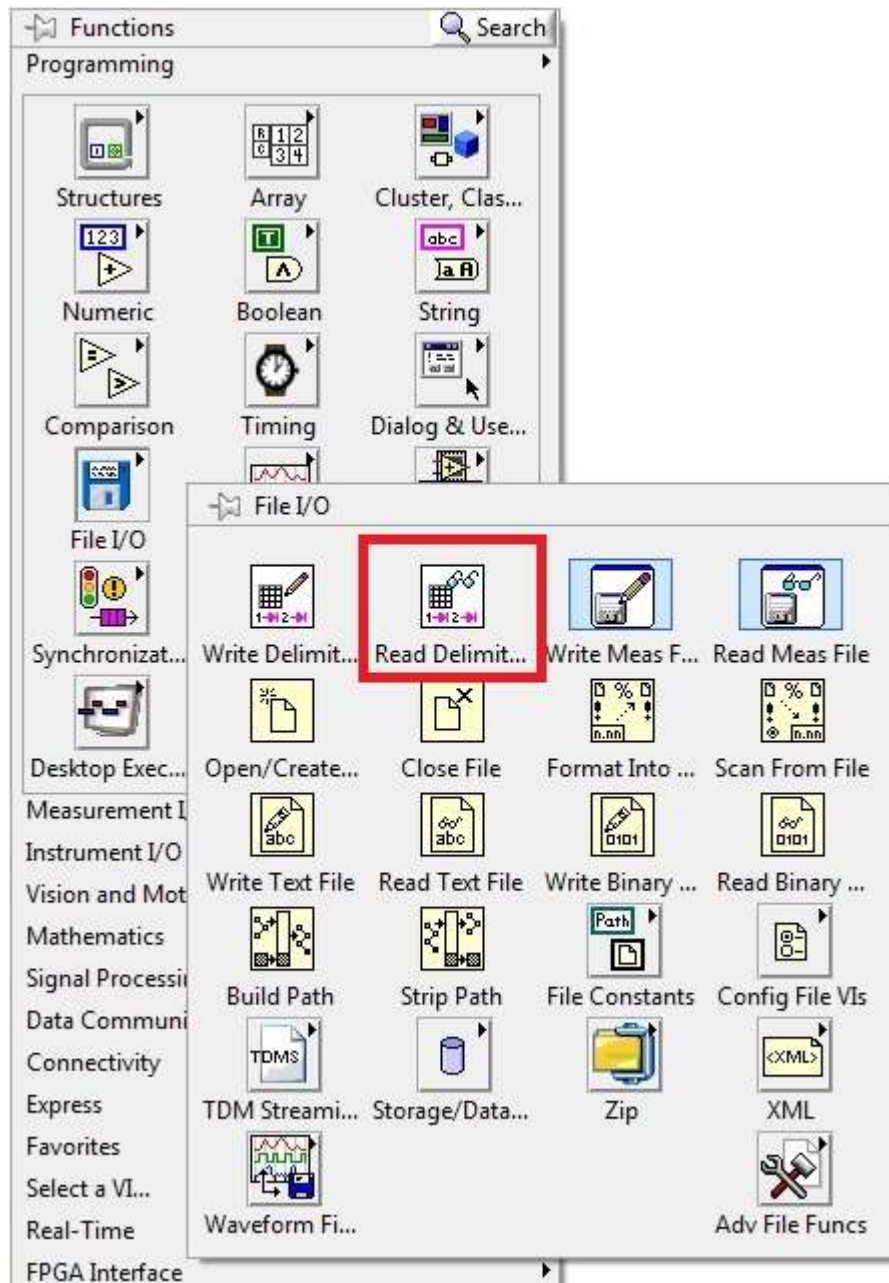

4. Right-click on the **delimiter (\t)** input located on the bottom edge of the VI and select **Create » Constant**. The default value for this input is a single tab character (`\t`).

1. For a comma-separated value file, change the newly created string constant to a comma.

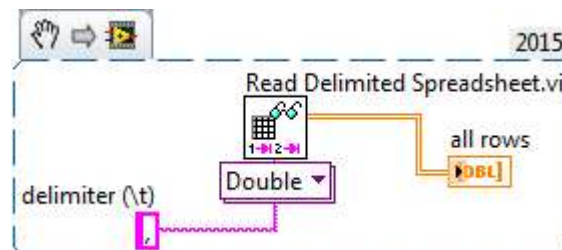

By default, the VI will display a File Dialog box to locate the `.csv` file at run-time. Additionally, you can display the data in multiple ways. The above example displays all of the rows read from Excel.

## Additional Information

**Note:** The Read Delimited Spreadsheet sub-VI contains a significant amount of functionality. Refer to the Context Help for more information about the various inputs and outputs.

To use the Excel file (instead of a CSV file) directly in LabVIEW, the Report Generation Toolkit is required. As of LabVIEW 2014, this toolkit is now included with LabVIEW Professional. Please see the product page for more information: [NI LabVIEW Report Generation Toolkit for Microsoft Office](#)

In LabVIEW 2014 SP1 and earlier, this functionality was encapsulated in the Read From Spreadsheet VI, which has been deprecated as of LabVIEW 2015.

---

## Related Links

- [Reading Microsoft Excel Numeric Data Into LabVIEW](#)
- [Report Generation Toolkit Compatibility with Microsoft Office and LabVIEW](#)

## Other Support Options

### Ask the NI Community

Collaborate with other users in our discussion forums

- [Search the NI Community for a solution](#)

### Request Support from an Engineer

A valid service agreement may be required, and support options vary by country.

- [Open a service request](#)
- [Purchase or renew support services](#)